
mamivac®

THERMO-PACK

DE
EN
FR
ES

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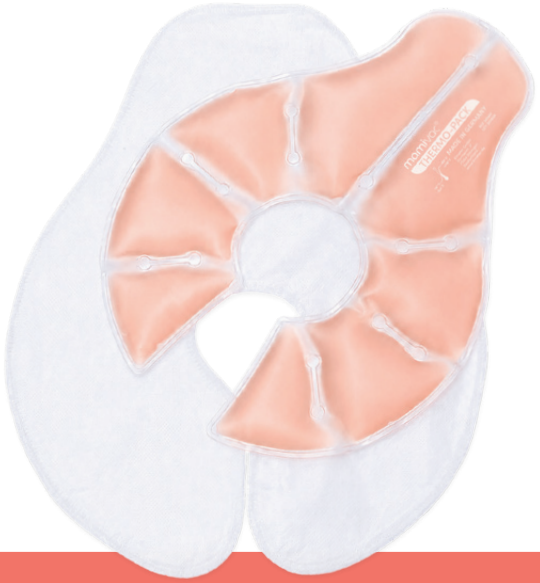
mamivac® THERMO-PACK

Instructions for use
Breast gel cushion

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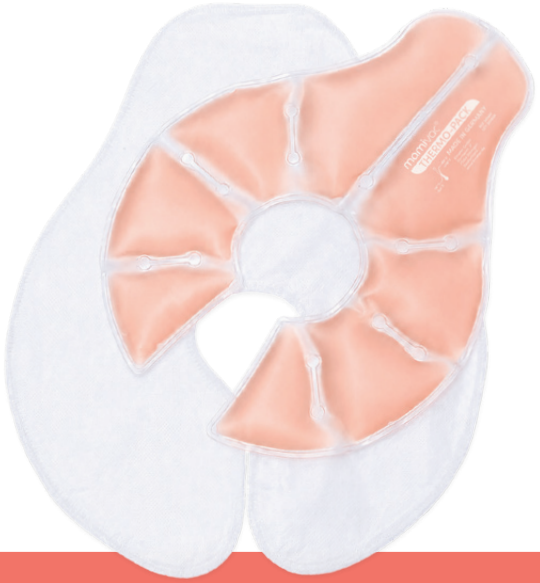REF 281487

Contents

1. Intended purpose 10
2. Design and components 11
3. Application 12
3.1 Correct positioning 12
3.2 Heat application 12
3.3 Cold application 12
4. Warnings and safety information 13
5. Cleaning and disposal 13
6. Product overview 13
7. Explanation of symbols 13

1. Intended purpose

The mamivac® Thermo-Pack is an anatomically shaped gel pad for warming and cooling the breast.

Before breastfeeding or expressing milk, warmth can stimulate milk flow. After breastfeeding or expressing milk, cooling can relieve feelings of tightness.

The mamivac® Thermo-Pack can also be used to warm or cool other parts of the body, such as the knee or neck.


2. Design and components

The mamivac® Thermo-Pack consists of an anatomically shaped gel pad with eight gel chambers and a matching protective nonwoven cover. The outer film of the gel pad is made of phthalate-free flexible PVC.

The gel chambers contain thermal gel that stores heat and cold. The central opening ensures that the nipple and areola remain uncovered when the pad is applied. The lateral extension reaches under the armpit.

The soft protective nonwoven cover is placed between the gel pad and the skin, protecting the skin from direct exposure to heat or cold. A hook-and-loop fastener keeps the cover closed.

a) Protective nonwoven cover (3 pieces)
b) Hook-and-loop fastener
c) Lateral extension
d) 8 gel chambers filled with thermal gel
e) Central opening

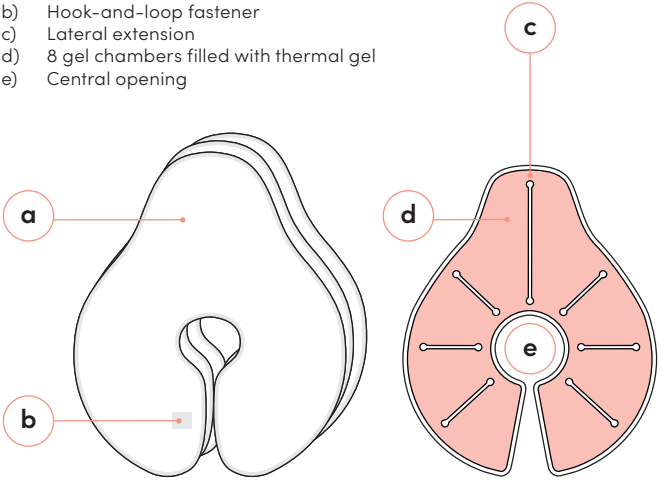Heating:
Place in a hot water bath (approx. 60 °C) for approx. 10 minutes.

Cooling:
Place in the refrigerator for approx. 1 hour.

Material:
The high-quality thermal gel withstands temperatures from –18 °C to +80 °C and is non-toxic.

3. Application

3.1 Correct positioning

Position the Thermo-Pack on the breast so that the nipple and the entire areola remain uncovered in the central opening.

The lateral extension should point towards the armpit.

⚠ **Before each use, check the temperature of the Thermo-Pack with the back of your hand or the inside of your forearm. Temperatures that are too high or too low can cause skin injuries.**

Only use the Thermo-Pack with the protective nonwoven cover to protect the skin from direct exposure to heat or cold.

3.2 Heat application

Warm the mamivac® Thermo-Pack in a warm water bath for about 10 minutes. A water temperature of approximately 60 °C is sufficient. Place the warmed Thermo-Pack in the protective nonwoven cover and apply it to the breast for approximately 5–10 minutes immediately before breastfeeding or expressing milk.

Do not expose the Thermo-Pack to temperatures above 80 °C.

Never heat the Thermo-Pack in a microwave.

Failure to follow these instructions may damage the gel pad or cause it to burst.

3.3 Cold application

Place the mamivac® Thermo-Pack in the refrigerator for about one hour. Alternatively, it can be cooled in the freezer at temperatures down to –18 °C. However, we recommend cooling it in the refrigerator to avoid excessive cold exposure to the skin.

Place the cooled Thermo-Pack in the protective nonwoven cover and apply it to the breast for approximately 10 minutes after breastfeeding or expressing milk.

4. Warnings and safety information

Stop using the mamivac® Thermo-Pack if the gel pad is damaged or leaking. The thermal gel it contains is non-toxic but must not be swallowed.

The Thermo-Pack is not intended for use on infants or children. Keep it out of the reach of children.

Do not apply the Thermo-Pack to open wounds.

5. Cleaning and disposal

Clean and disinfect the mamivac® Thermo-Pack only with products suitable for flexible PVC.

Dispose of damaged or leaking Thermo-Packs in accordance with local household waste regulations.

6. Product overview

Thermo-Pack including 3 protective bags	REF 281487
Protective bags, 25 pieces per bag	REF 281498

7. Explanation of symbols

	Follow the instructions for use		Manufacturer
	Mandatory action sign		Warning sign
	Prohibition sign		Maximum temperature: +80 °C Minimum temperature: –18 °C

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