mamivac® EASY

DE
EN
FR
ES
AR

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mamivac® EASY

Instructions for use
Manual breast pump

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REF 281360

Before using the mamivac® EASY manual breast pump, read the instructions for use carefully and keep them for future reference.

Content


1. Notes on the instructions for use 12
2. Assembling and components 13
3. Purpose 14
3.1 Indication 14
3.2 Contraindication 14
3.3 Side effects 14
3.4 Target group of patients/users 14
4. General notes 14
4.1 Functioning of the manual breast pump 15
5. Hygiene and cleaning 15
6. Practical tips 16
6.1 Storing breast milk 16
6.2 Accessories 16
7. Classification 17
8. Disposal 17
9. Spare parts 17
10. Warranty 17
11. Technical data 17
12. Legend 17

1. Notes on the instructions for use

These instructions for use are valid for the mamivac® EASY manual breast pump as well as all other derivates of identical design.

mamivac® EASY	REF 281360
mamivac® EASY in poly bag	REF 281360B

2. Assembling and components

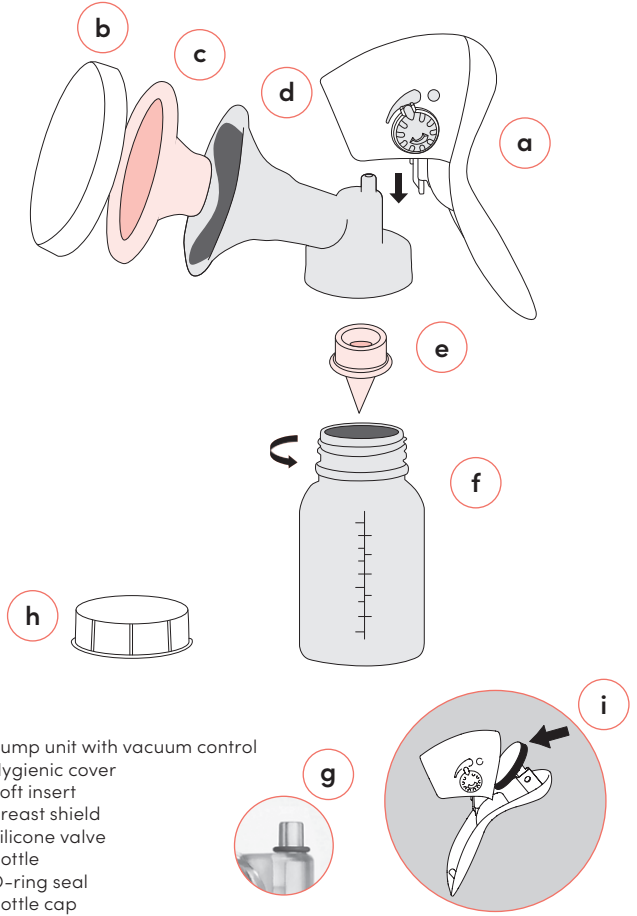a) Pump unit with vacuum control
b) Hygienic cover
c) Soft insert
d) Breast shield
e) Silicone valve
f) Bottle
g) O-ring seal
h) Bottle cap
i) Piston seal ring

3. Purpose

The manual breast pump mamivac® EASY is used in order to genty and softly express and collect breast milk.

3.1 Indication

Breastfeeding problems that require the use of a breast pump may be: mastitis (inflammation of mammary glands), aching or sore mamillae (nipples), anatomic variants of the mamillae which make feeding the infant difficult (flat or inverted nipples), insufficient or excessive milk production, premature delivery, metabolic disorders of the infant, anatomic disorder in the face of the infant, poorly sucking infant, temporary illness of the mother (medication intake).

3.2 Contraindication

Taking medication in general may affect breastfeeding an infant. The patient information leaflets of the pharmaceutical companies contain appropriate information.

mamivac® EASY may not be used:

- to suction other body fluids, such as blood, secretions or mucus
- to suction outside the medical field
- on other living beings.

3.3 Side effects

In individual cases, the skin may be irritated resulting in redness or red spots. In rare cases, sore nipples, bruising, or injuries to the nipples may occur.

3.4 Target group of patients/users

Breastfeeding mothers

4. General notes

Check the product for visible damage and proper functioning before each use. The pump must not be used further if visible damage (e.g. cracks, material fatigue, deformation, leakage) or non-removable residues are detected.

- Wash your hands thoroughly with soap before each expression and wipe the breast with a damp cloth. Dryness and cracking of the nipples can be avoided by not using cleaning agents such as soap or alcohol.
- Prepare yourself calmly for expressing. Sit comfortably and relax.


Hold the breast shield centrally against the breast with slight pressure so that it seals properly. Begin expressing.

- The pump is equipped with a stimulation mode that promotes milk flow. To support this, you can additionally massage the breast with your other hand
- Express from both breasts alternately, unless instructed otherwise by the doctor. The duration of expression and the amount of milk obtained vary idividually and depend on many factors

- If you experience pain while expressing or can express only a small amount of milk, please consult your midwife/lactation consultant or another healthcare professional.
- The use of the breast pump should be avoided in the presence of open wounds or severely damaged nipples.
- All serious incidents related to the product must be reported to the manufacturer and the relevant national authority.

4.1 Functioning of the manual breast pump

The stimulation phase (1):
To trigger the let-down reflex, turn the vacuum control to the very left.
Press the pump handle dynamically downwards several times.
This imitates the baby's natural sucking behaviour. Keep the pump handle pressed for 3-4 seconds, so that the milk starts to flow into the bottle.

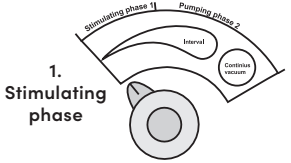The expressing phase (2):
When the milk is flowing, turn the vacuum control to the right. You can express now by pressing the pump handle downwards and keeping the handle pressed until the vacuum has been released. If the milk does not start to flow right away, relax and try again – even at a later point in time. If expressing causes pain, do not use the pump anymore! Please contact your doctor or your midwife/lactation consultant. Make sure not to express too much milk into the bottle. Keep the pump in an upright position.

5. Hygiene and cleaning

Before you use the breast pump for the first time, clean components b) to h) and boil in boiling water for at least 5 min. or use a steam disinfection device, such as a steriliser for instance. Please observe the manufacturer's corresponding notes.

Do not use a microwave oven. 

- Do not use chemical disinfection baths, solvents or abrasives.
- All components b) to h) can be cleaned with commercial washing-up liquids by hand or in the dish washer.
- Prior to using the pump the next time, clean the parts which came into contact with milk, and check them for proper functioning and condition.
- After cleaning, rinse components b) to h) thoroughly with lukewarm water to remove residues of the washing-up liquid, especially on those components which come into direct contact with breast milk. You might use a brush to clean areas which are difficult to reach.
- Assemble the manual breast pump according to the drawing in »Assembling and components«.

6. Practical tips


Wrap the cleaned and dry parts in a clean cloth and store them in a container when not in use.

6.1 Storing breast milk

Breast milk can be stored as follows:

- Freshly expressed at room temperature (16-29°C) 3-4 hours
- In the rear of the fridge (below 4°C) 3 days
- In the freezer (-18°C) 6 months

Fill the used bottles or containers to two thirds only and write down the fill date. Instead of a bottle, you can also use our specially designed mamivac® breast milk freezer bags [REF 282100]. Do not defrost frozen breast milk in the microwave, but at room temperature. Place the bottle or bag in warm water (37°C) as required. After defrosting, gently shake the container or bag to mix the milk.

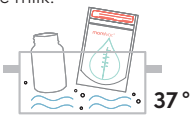6.2 Accessories

Soft insert	REF 281215
mamivac® breast milk freezer bags PU = 20 pcs.	REF 282100
Hygienic cover	REF 49010147

7. Classification

Medical device of risk class I

8. Disposal

Please dispose of your mamivac® EASY manual breast pump in accordance with your local disposal regulations.

9. Spare parts

Silicone valve:	REF B49054105	Hygienic cover:	REF 49010147
Piston seal ring:	REF 49054125	Standard breast shield:	REF 28012011
O-ring seal:	REF 49054005	Bottle 120 ml with cap:	REF 28011131
Retaining spring:	REF 49040100		

10. Warranty

Ignoring the notes in these operating instructions or an improper use of the device will void any warranty claims. Kimetec GmbH rejects any liability for personal injuries and material damages arising from this.

11. Technical data

Max. vacuum: -33 kPa (±15 %)
H x W x D: 166 x 152 x 75 mm
Breast shield: Ø 26 mm
Weight: ca. 110 g
Soft Insert: Ø 24 mm

12. Legend

	Consult instructions for use		Manufacturer
	Safety notice		Warning
	Prohibition sign		Keep away from sunlight
UDI	Unique Device Identifier	REF	Catalogue number
LOT	Batch code	MD	Medical device
	Transport and storage: -25°C to +70°C \| Operating conditions: +5°C to +40°C		
CE	CE marking according to Regulation (EU) 2017/745		

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